

F-300DT

Digital Facial Steamer

User Manual

Welcome to use this product! A0
Please read this instruction carefully before use.

About This Manual

This manual shall only be used for the reference of operation of this particular model. If the manual is used for other purposes, we will not undertake any according consequences and responsibility. All information of this manual shall be protected by copyright and shall not be copied and translated in any part of this specification without the written consent of manufacturer.

Safety Guidance

To ensure that you use the product correctly, please read all the security guidance carefully before use. This security guidance can help you use this product correctly to avoid harm to you and your family, and unnecessary property damage. When using this product, especially when a child is present, basic protective measures must be taken. The security guidance are as follows.

	Warning!	Indicates a dangerous situation that may occur. If not avoided in time, it may lead to serious injury or even death. • This product must be used in accordance with the inherent purpose stated in this manual and cannot use the accessories, which are not recommended by manufacturer. • The following people should not use: ① Women during pregnancy or menstruation, ② People with allergies, skin diseases, or severe sunburn, ③ People with skin damage, eczema, swelling, etc. ④ People with limbs, senses, or mental illness, ⑤ People who have just completed surgery and are recovering. • When this product is working, do not wind the power cord around the body, and do not use an extended power cord. Keep out of reach of children. Do not allow children to use this product. • When using and placing, the power cord must be kept away from high temperature. • When this product is powered on, it must not be left unattended. • When the power cord and plug of this product are damaged, work abnormally, overturned or dropped in water, you cannot continue to use it. You should send it to the service center for inspection and repair. • This product is not recommended for children. • Do not insert or put any foreign objects into the body and nozzle of this product.
	Attention!	Indicates a dangerous situation that may occur. If it is not avoided in time, it may cause personal injury and property damage in severe cases. • During non-use, if you need to move the steamer with hot water, be sure to handle it carefully. • Please keep this product properly and keep it out of reach of children. • When using, do not remove the water cup. • Do not touch the herbal cap during work.
	Dangerous (avoid electric shock)	• Unplug the power cord immediately after use. • Do not use this product while bathing. • Do not place or store this product in a location that can easily fall into the bathtub, washbasin, sink, etc. • Do not allow water or other liquids to enter the product due to improper operation. • If the product falls into the water accidentally, please do not touch the product and unplug the power cord immediately. • Do not use this product in a humid environment.
	High temperature attention	When working with this product, please pay attention to high temperature parts such as nozzles and funnel cap to avoid burns.

Product Introduction

This product is used for supporting other beauty instruments and beautician trolley.
The models are F-300D / F-300H / F-300E / F-300F / 3021 / 3051 / 3052.

F-300H



F-300E




3051




3021



Components

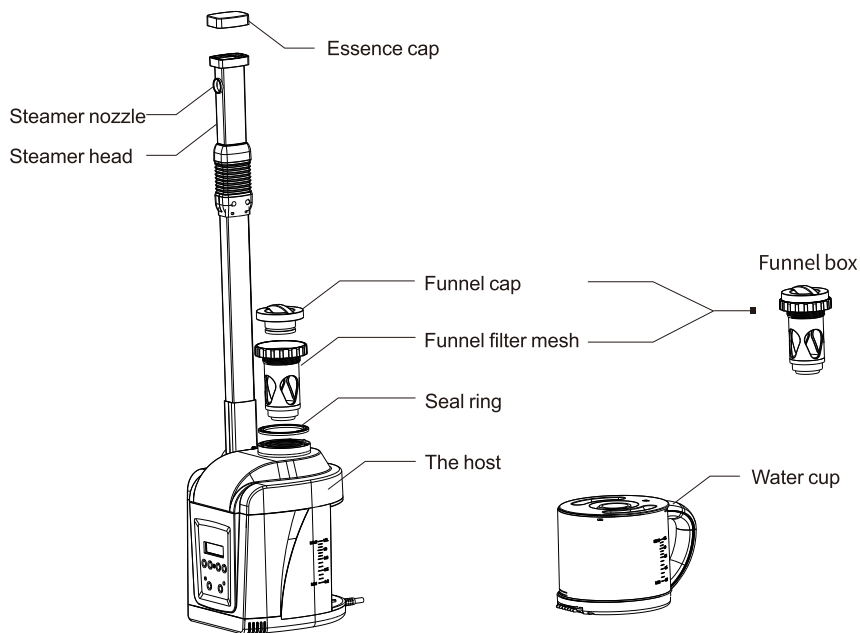

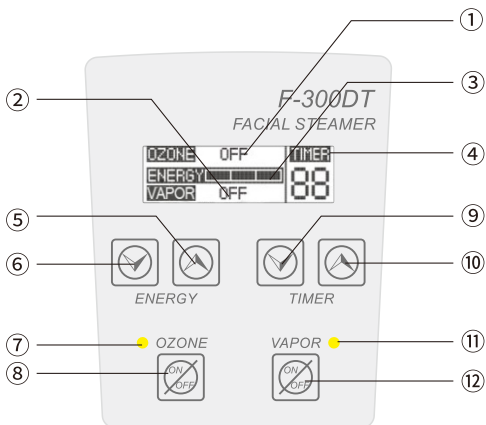

- ① Ozone display ② Spray display
- ③ Energy ④ Timer
- ⑤ Energy increasing button
- ⑥ Energy decreasing button
- ⑦ Ozone indicator light
- ⑧ Ozone switch ON/OFF
- ⑨ Time decreasing button
- ⑩ Time increasing button
- ⑪ Vapor indicator light
- ⑫ Spray switch ON/OFF



When you use this product for the first time, you may find that there is residual water around the water cup. This is the distilled water remaining in the factory used in the aging test, which will not cause any hygiene problems.

Product Parameters

Model Type	F-300DT	
Frequency	50/60Hz	
Wattage	650W/800W	
Voltage	AC 100-120V	AC 220-240V

Operation Method


● Add water Method One

First, unscrew the funnel cap(picture ①). Then press and hold the gray button on the host(picture ②), and take out the water cup. Add water to the water cup. After adding water, it needs to be pushed into the pit of the bayonet when it is inserted into the host. After it is inserted, check whether the water cup is installed in place.


● Add water Method Two

First, remove the funnel cap and add water. After adding water, please make sure that the funnel cap and funnel filter mesh are tightened.

Please make sure that the water in the water cup is between the safe water level (MIN-MAX).

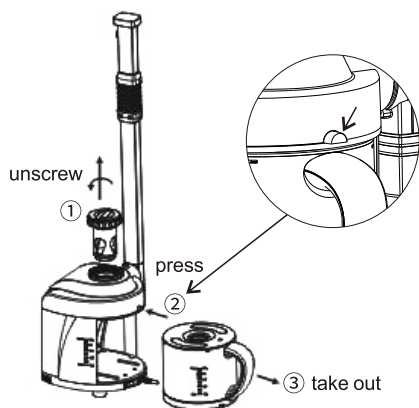

Add water Method One

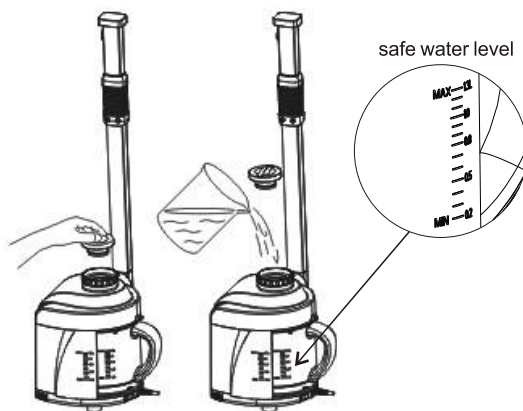

Add water Method Two

- Connect power supply
- After the host is powered on and self-checked, it enters the standby state.
- Adjust time and energy.
(When the water is at room temperature, the energy should be adjusted to the maximum to shorten the heating time) .
- After the water is boiled, press the ozone switch then the ozone indicator light is on. The steam is obviously converted into a fine mist and sprayed out with smell of ozone. At this time, it enters the normal spray state and the beauty operation can be performed.
- For the function of aromatherapy spray, please open the essence cap on the top of the nozzle and add 2 drops essence. Then close the essence cap (as picture 1)



Picture 1



• During operation, do not move water cup and do not tilt, shake or move this product. Otherwise hot water may splash out or a large amount of the steam is ejected from the nozzle.



• When working with this product, please pay attention to high temperature parts such as nozzles and funnel cap to avoid burns.

Cautions

Operation Cautions



Warning


- Do not start up when the water level is lower than the lowest water mark.
- Do not use this product under sun for long term. Otherwise it may cause unstable performance.
- Do not disassemble the shell, host or other parts without permission. In case of abnormal work, please contact the dealer for professional guidance. Or ask professional maintenance personnel for inspection and maintenance.




Attention

- When the water in the water tank is below the lowest water level, the host will emit a beep tone and stop spray. At this time, should add proper water and restart the heating function to work normally.
- Can use purified water or tap water. Do not use distilled water. Do not add viscous substances such as flavors or herbal to the water cup. Otherwise, it is easy to cause water spray and cause damage to the machine. When the machine stops working, please turn off the power.
- The ozone has a bactericidal effect. But long-term use will cause damage to the skin and destroy the nutrition of the herbal. The use time should be selected according to different skin types.
- During operation, do not spray directly on the skin. It may cause burns. You must wait until the nozzle can spray the fine mist stably and evenly before operation.
- If there is water leakage during work, first take out the water cup to check whether the seal ring is well set and then reinstall the water cup. If the situation still does not improve, please contact the dealer to resolve.


- **When using this product for the first time, please use tap water, and then use pure water, which can try to avoid the generation of the scale.**
If the detection function responds incorrectly, it may be that the water quality is too pure. At this time, it can be solved by working with tap water once.
- **This instrument has a detection circuit. If there is no water in the water cup, the funnel cap is not tightened, etc., the host may not work. At the same time, the host will also emit "beep" "beep" prompts. The above problems should be eliminated before starting the host.**

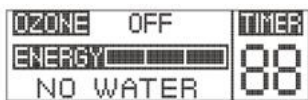

No water

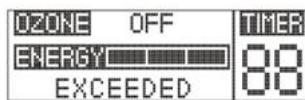

Overfilled

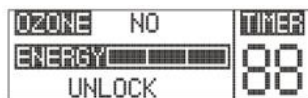

The funnel filter mesh is untightened

- Please keep this facial steamer away from the high-interferential device (like: high-frequency instrument). If the circuit board or LCD screen is accidentally interfered, the steamer may stop working suddenly or the quality of performance will be lowered. Normally, please turn off the power, then re-start after 15 seconds. This product can be recovered to the normal working status. If still not improved, please unplug the power immediately and contact the dealer for repair.
- If the main body of the funnel cap is not tightened, it will cause steam leakage. But too tight will make the sealing ring of the atomization chamber easy to age, resulting in air leakage. To avoid air leakage caused by the aging of the seal ring of the vapor dispenser, please change the seal ring in every 3 months. Each product has one more seal ring as spare. If it is lost or damaged, please contact the dealer to get a new one for replacement.
- Do not add water without the funnel filter mesh. It may easily allow the water to flow to other parts and cause a short circuit.



Do not add water without the funnel filter mesh

- Do not pull out the water cup without removing the funnel filter mesh. Otherwise, the funnel filter mesh will be broken.
- Do not remove the host by grabbing the spray arm to avoid loosening and damage to the spray arm.
- Please note, take out the water cup after using it every day, and check whether there is any water drops or mist in the place where the water cup is connected to the power supply. If it is, please clean it immediately with a dry towel and keep the host dry and ventilated. Please turn off the power before cleaning to avoid damage to the product or electric shock.
- During operation, you must wait until the fine mist is sprayed stably and evenly before it can be directed at the area that needs to be sprayed. The corresponding distance must be maintained when spraying. The distance can be compared with the following reference table.



Do not pull out the water cup without removing the funnel filter mesh

Different skin types refer to different distances. The following data is for reference only.

Spray time reference table		
Skin Type	Distance	Time
Oil skin, poor blood circulation skin	35-40cm	15min
Normal skin	35-40cm	15min
Sensitive and dry skin	40-45cm	5min
Inflamed skin	40cm	5min

Instrument Maintenance

- When operating, please pay attention to the water level of the water cup.
- After operation, please clean the product in time.
- The scale on the water cup and heating plate needs to be cleaned at least once a week or after using funnel twice. Can be scrubbed with a soft metal brush or soaked in vinegar for 5-7 minutes. Then rinse it off with tap water. Please note that funnel can not be reused.
- Check the seal ring every 2 weeks. If it is deformed or deteriorated, please replace it in time.

Troubleshooting

- When an abnormality occurs, stop using it immediately and cut off the power supply.
- After normal power-on, if the product heats up for no reason before the spray function is turned on. Please cut off the power supply and contact dealer for repair.

Abnormal phenomena	Cause	Solution
The sprayer does not work when press the switch button	1. Power is not connected	Connect the power supply
	2. Do not adjust the energy gear	Press energy increasing button
The sprayer does not work	1. Not powered on	Boot according to the operation method
	2. Water shortage	Add water according to the operation method
	3. Nozzle clogged	Clean the product according to the maintenance guideline

If problem is still not solved after above operations, please contact maintenance center or supplier to solve it.



Do not dispose of old or defective appliances in domestic garbage; this should only be done through public collection points.

- ▲ Please check the product specification before use. This manual is for reference only. Specifications are subject to change without prior notice.